




Lavender
Medical

IO FiX™

Fusion Fixation

The IO FiX™ Advantage: *Designed for Fusion*

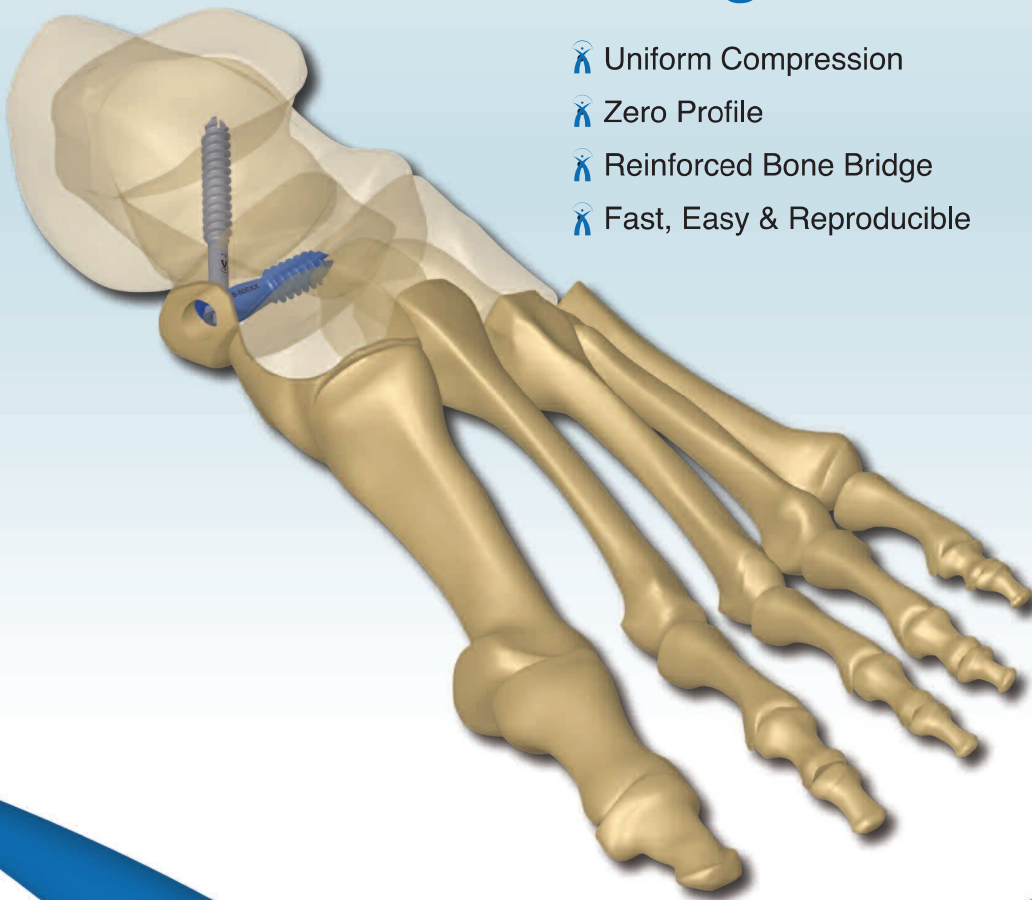

- Uniform Compression
- Zero Profile
- Reinforced Bone Bridge
- Fast, Easy & Reproducible

Because Uniform Compression
is **Better** Compression


EXTREMITY
MEDICAL™

How IO FiX™ Works

Zero Profile

Implants are placed within the bone



- ⌘ Minimizes soft tissue irritation
- ⌘ Decreases need for hardware removal

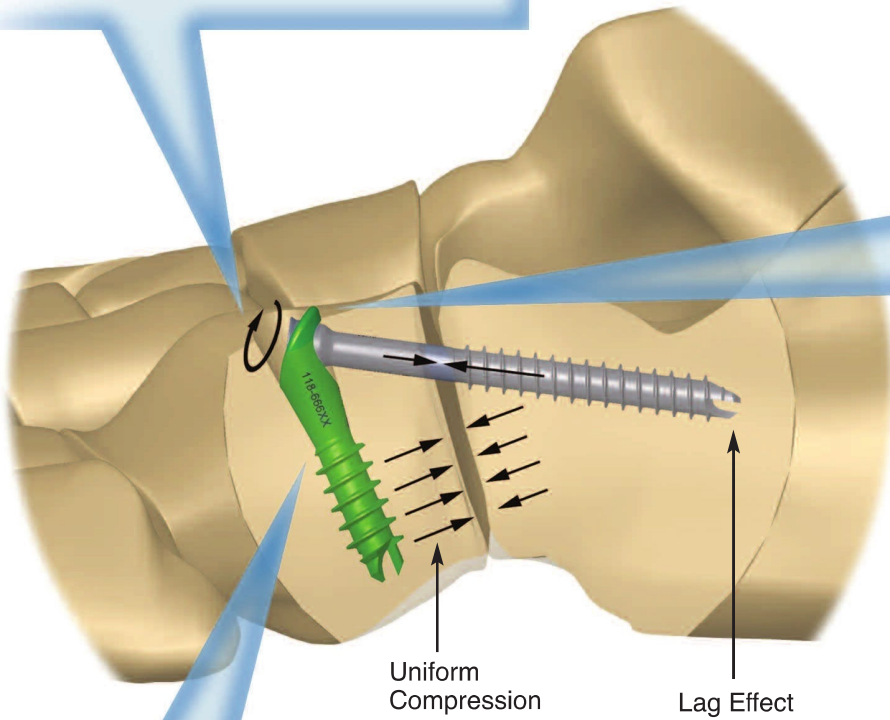

Reinforced Bone Bridge

Screw lags against a reinforced metallic bone bridge

- ⌘ Unlike screws, IO FiX™ maintains compression if the cortical bridge is compromised

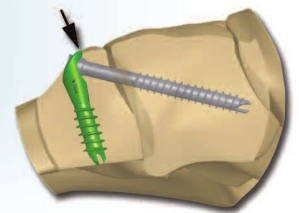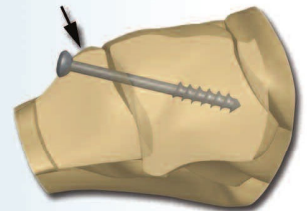


X-Post™

Distributes compressive forces across a greater surface area

- ⌘ Unique Morse Taper Locking System lags and locks simultaneously
- ⌘ More uniform compression*
- ⌘ Higher peak compression*

Designed for Fusion


Uniform Compression*

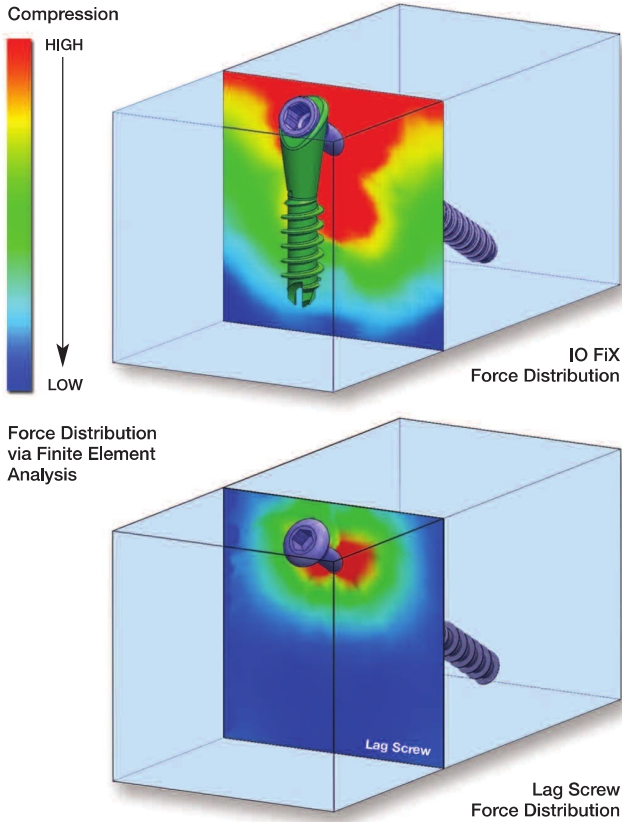


Stability**

Superior fatigue and strength resistance decrease the likelihood of a construct failure


Fatigue Strength*

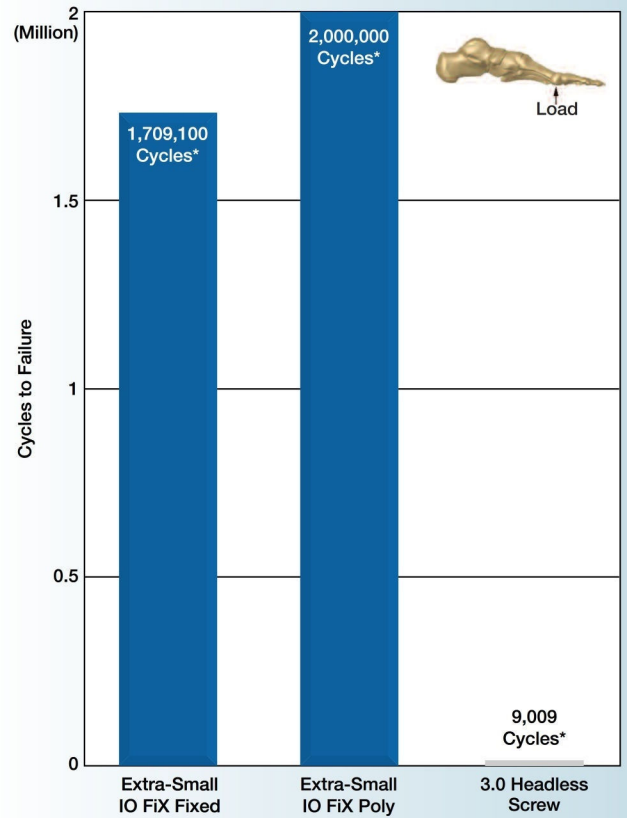

Peak Compression*

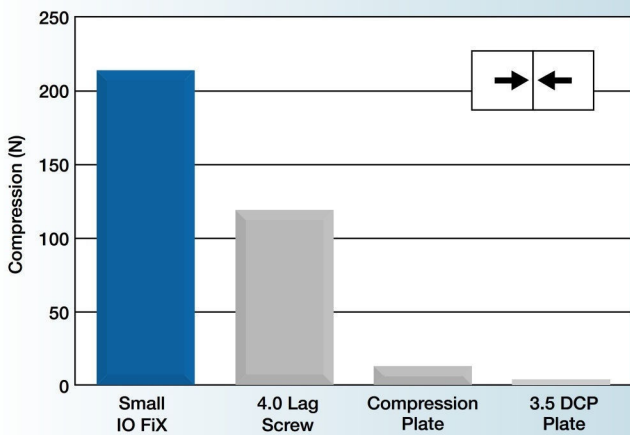

Bending Strength*

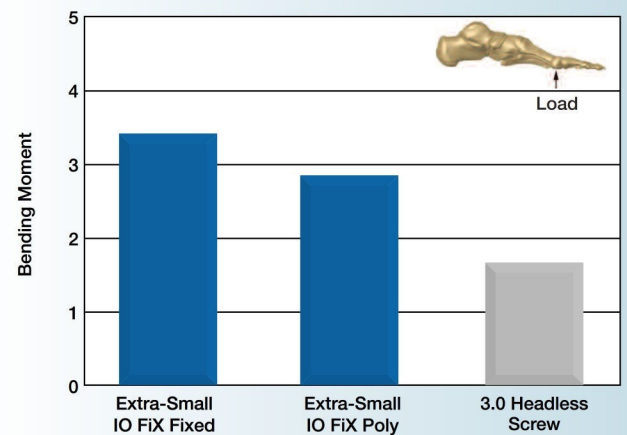

* Data on file, Extremity Medical ** Manual of Internal Fixation - A

Versatility

Various sizes available for many foot and ankle indications



Lapidus



Talonavicular Arthrodesis



Triple Arthrodesis



Ankle Fusion



MTP Fusion



TMT Fusion



Calcaneal Osteotomy

300 Interpace Parkway
Suite 410
Parsippany, NJ 07054
Phone: 973.588.8980
Customer Service: 888.499.0079
Fax: 888.499.0542
www.extremitymedical.com



Distributed by:
Lavender Medical Limited
tel : 0845 676 9733
email : info@lavendermedical.com
web : www.lavendermedical.com



Patent Pending
Extremity Medical and IO FIX are trademarks
of Extremity Medical, LLC.
©Extremity Medical, LLC. 2012. All rights reserved.
LBL 118-99105-EN Rev B